







Silensor-sl

Verbinderswechsel
exchange of a connector



45 °C
113 °F






• changement d'attache • cambio de conector • sostituzione congiuntore • substituição do conector

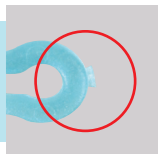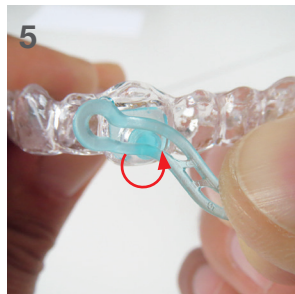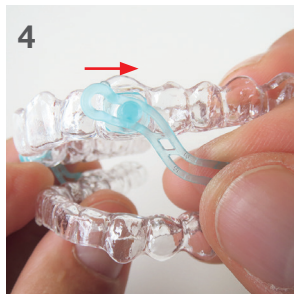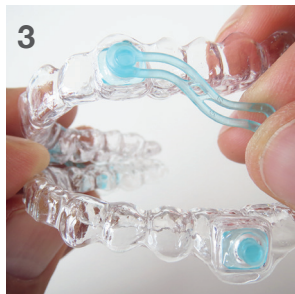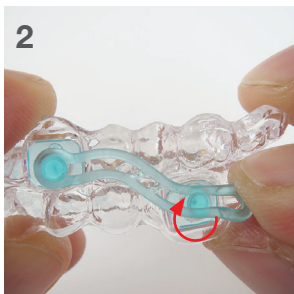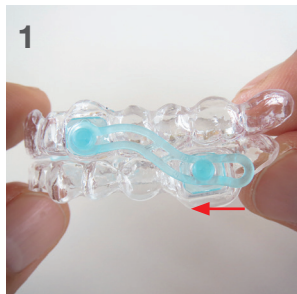

- loser Anker
- loose anchor
- ancre lâche
- ancla suelta
- ancoraggio staccato
- soltar suporte





- Anker halten
- hold anchor
- tenir ancre
- mantener ancla
- tenere premuto l'ancoraggio
- segurar suporte

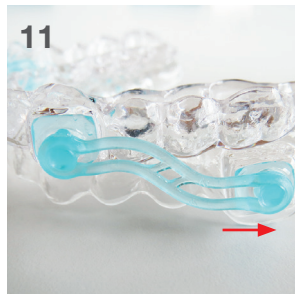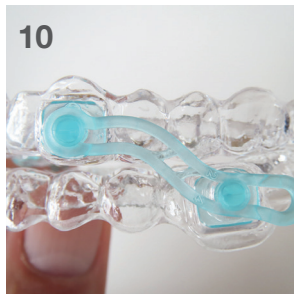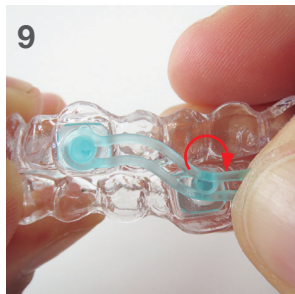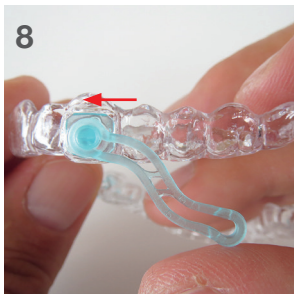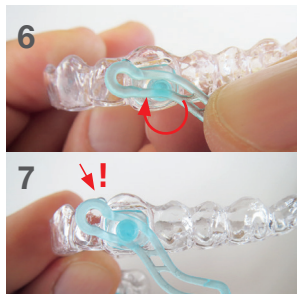


restliche Verbinder hier aufkleben
glue remaining connectors here
coller attaches restantes ici
pegar conectores restantes aquí
colla premuto l'ancoraggio qui
colar aqui os restantes conectores



zuerst verwendeter Verbinder
first used connector
attache en premier utilisée
conector utilizó por primera vez
congiuntore prima utilizzazione
conector a utilizar na primeira vez

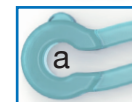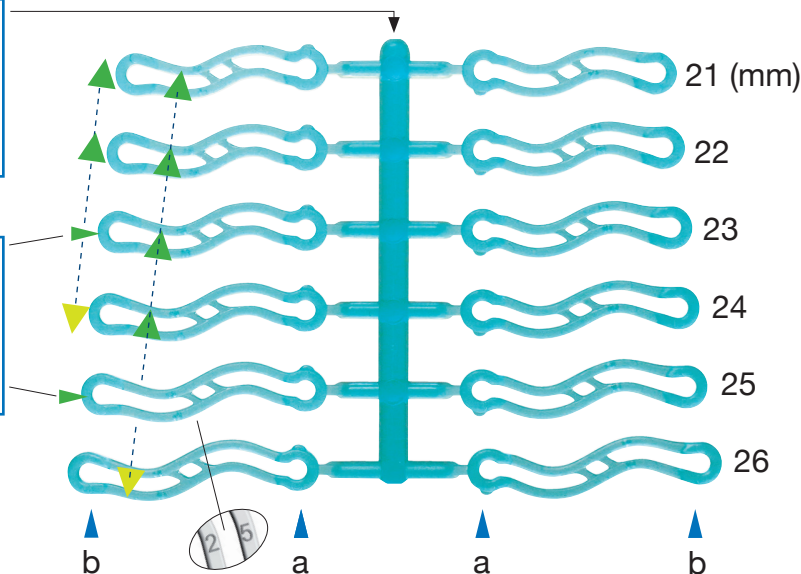

oberer Eckzahn
upper canine
canine supérieure
canino superior
canino superiore
canino superior

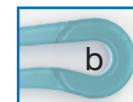

unterer Molarenbereich
lower molar area
région molaire inférieure
zona de molares inferior
zona dei molari inferiori
zona do molar inferior

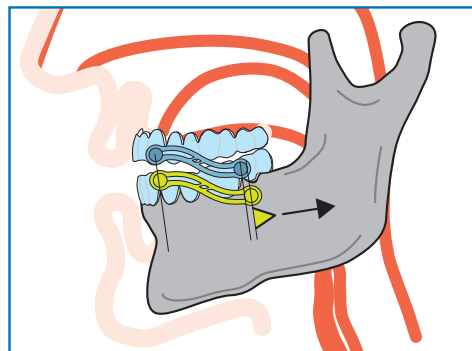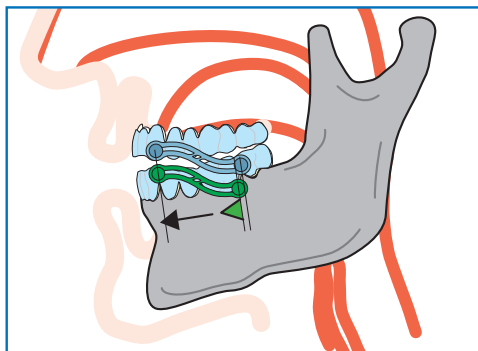

sl-Zange/sl-pliers/sl-pince/sl-pinza
Empfehlung zum leichteren Verbinderwechsel/
Recommendation for an easier connector
replacement, Ref.: 59 60 60



DE



EN



FR



IT