


25 colours are available:



Please indicate colour numbers: ☒

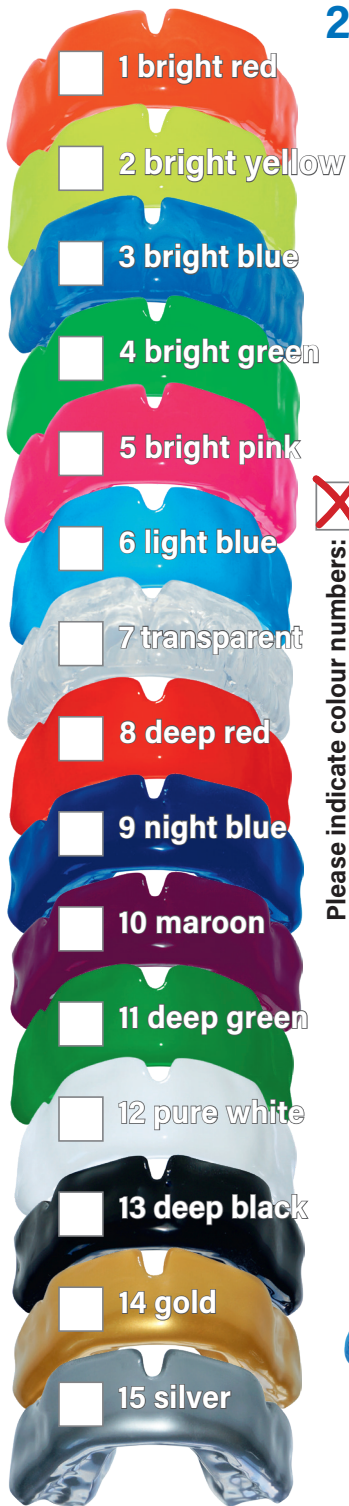

color deviations possible

Playsafe 4u

Free individualisation of your Playsafe triple




Playsafe triple mouthguard out of up to **4** colours!

You can compose your own multi-coloured mouthguard out of the colours **1** to **15**.

Please indicate the colour numbers:

Example: →

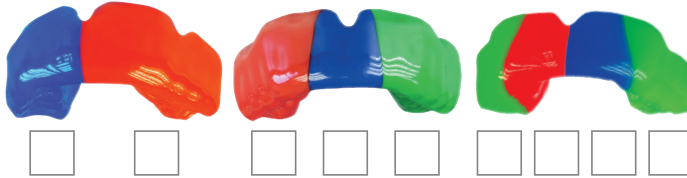

Please check the requested Playsafe triple type: ☒


PLAYSAFE® triple + **PLAYSAFE® triple-light**



Playsafe triple and triple light are suited for any kind of sport.

Playsafe triple light is suited in case of cramped oral situation.

Recommended by 'Verband der Deutschen Boxringärzte e. V.'

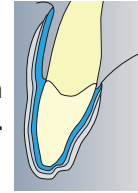

Dear Athlete,

Every Playsafe triple mouthguard is fabricated for you by a trained person according to the Playsafe guidelines.

You will get your Playsafe mouthguard in your dental practice. There an impression of your upper and lower jaw will be taken.

In the dental lab, the plaster model and the sports mouthguard will then be produced.

The control of fit and handing over will again take place in the dental practice.

An original Playsafe triple sports mouthguard is always 3-layered and shows the Playsafe triple type and the name of the athlete.

Your Playsafe triple team:

ERKODENT®

Erkodent Erich Kopp GmbH • Siemensstraße 3 • 72285 Pfalzgrafenweiler
Germany • Tel. + 49 (0) 74 45 / 85 01-0
info@erkodent.com • www.erkodent.com
EN ISO 13485 • ISO 9001

PLAYSAFE® triple

Sports mouthguard



CE The first individually fabricated sports mouthguard with passed EU type examination



ERKODENT®




Why individually fabricated Playsafe triple sports mouthguards?

Playsafe triple sports mouthguards fit exactly onto the teeth, do not move and allow to breathe freely and to speak normally.

Sports mouthguards are not all equal!

Boil & Bite mouthguards fit poorly and become loose. The athlete has to hold this mouthguard onto the teeth, this prevents him from producing outstanding performances.

Many scientific studies show that a Boil & Bite mouthguard provides a false sense of protection due to the dramatic decrease in thickness when the athlete bites it into place during its softened state (biting through).

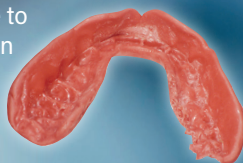

Boil & Bite mouthguard

First custom-made sports mouthguard with **CE** marking to denote their compliance with the essential health and safety requirements of Regulation (EU) 2016/425 as 2nd Category PPE = Personal Protective Equipment against medium risks.

3-layered
PLAYSAFE® triple-light
Sports mouthguard, flexible with hard inter-layer:


▶ hard ▶ soft, flexible

Playsafe triple (-light) is suited for all sports.

The hard inter layer comprises the whole guard.

Highest protection with smallest extent!

(total thickness Playsafe triple 5.5 mm, Playsafe triple light 4.1 mm)

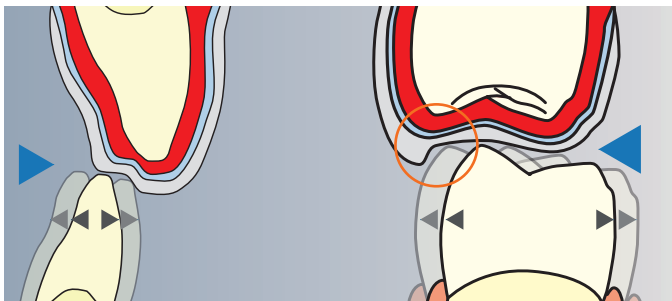

Playsafe triple (-light) shows a „buccal shield“ ▶ despite the movability of the antagonist.

This circulating protective border (orange circle) prevents too large movements of the lower jaw in case of lateral impacts.

Tests: intensive endurance tests show the high protective effect of Playsafe triple mouthguard types.

Resistance against deforming or pendulum test to measure power absorption and distribution.



... by the way, every Playsafe triple mouthguard can be adjusted to changing proportions of teeth.

https://www.erkodent.de/wp-content/documents/products/PSRefix_EN_web.pdf

